

## Image capture and mark-up software

- View, capture, mark-up and document images – quickly and simply
- Easy control of camera settings
- Record .AVI video
- Multi-language software



FM 557119

Vision Engineering Ltd has been certified for the quality management system ISO 9001:2008.



## Image capture and mark-up software

Quickly and simply capture and document images with ViCapture software.

Designed specifically for use with compatible Vision Engineering products\*, ViCapture allows users to quickly optimise camera settings, capture an image and add any required mark-up.



All in a highly intuitive, easy to use software package.

### Highlights

- ✓ View, capture, mark-up and document images, quickly and simply
- ✓ Mark-up with text, arrows, shapes etc.
- ✓ Optimise and save camera settings for easy recall

### \*Compatible with...

- ✓ Mantis Elite-Cam HD eyepiece-less stereo microscope
- ✓ Lynx EVO eyepiece-less stereo microscope
  - with Smart Cam integrated HD camera option
  - with Smart Cam 5 integrated 5MP camera option

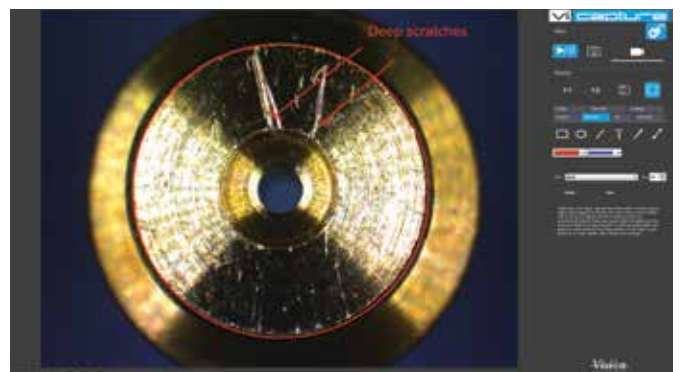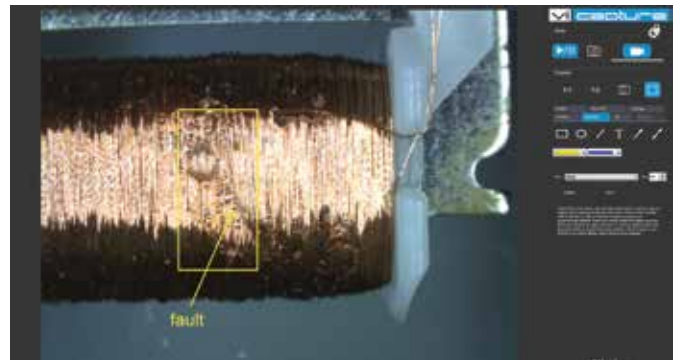

[Watch the video >>](#)



[www.visioneng.com/vicapture](http://www.visioneng.com/vicapture)

## Image capture & video recording, made easy

With large icons and tabbed menu panel, the software provides highly intuitive operation – even for first-time users.


Quickly optimise camera settings to view the best image, then simply click the image capture or video record icons. Settings can be saved for easy recall – for example, quickly recall optimised camera settings for highly reflective, light or dark-coloured subjects.

## Image mark-up tool

Important features can be clearly highlighted with a choice of shapes, arrows and text – in a selection of colours for best contrast.

## Features

- Easy control of camera settings
- Save images as JPG, BMP, or PNG files
- Mark-up with text, arrows and shapes
- Record .AVI video for up to 30 seconds
- Focus aid – displays a numerical value for peak focus
- Save optimised camera settings for easy recall (e.g. highly reflective subjects etc.)

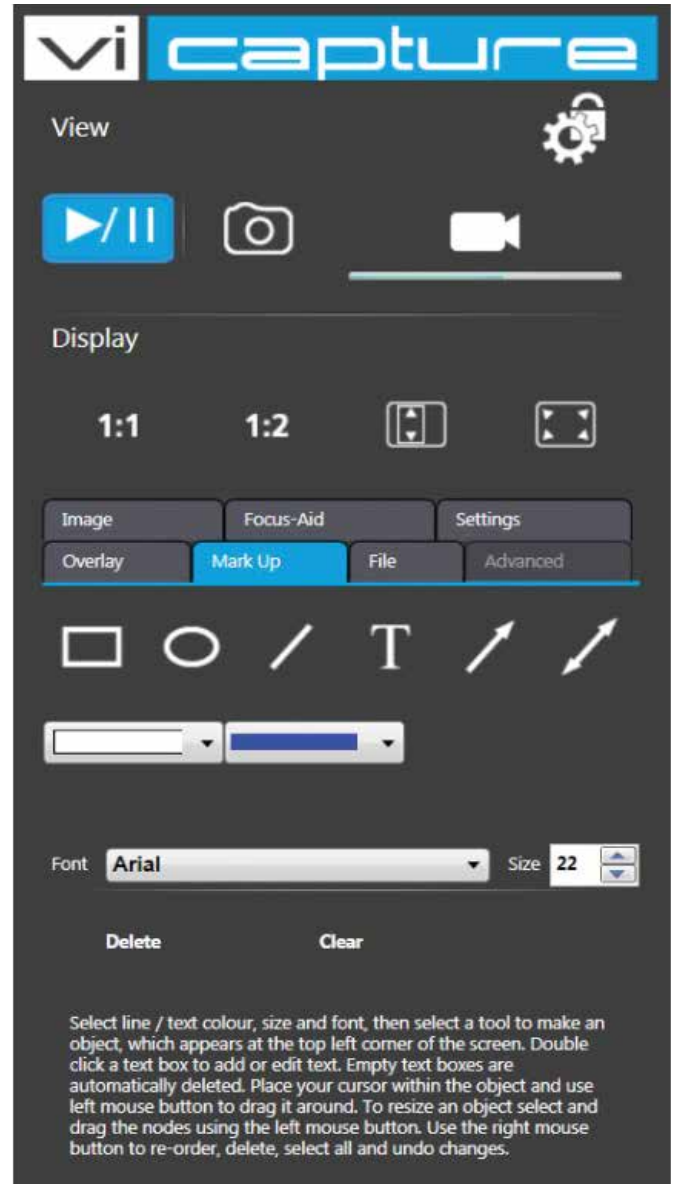

## Multi-language software

Supported languages: English, French, German, Italian, Spanish, Portuguese, Polish, Turkish, Russian, Chinese, Japanese, Korean

## System requirements

Windows 7 (or later), min. 500MB of storage, Intel i3 or better, 2GB of RAM, .NET framework and is compatible with Vision Engineering products\*

## Do you require on-screen measurement and reporting?

Vision Engineering has a range of software solutions offering enhanced functionality, including on-screen image measurement and advanced reporting.

**Find out more >>**

[www.visioneng.com/products/software](http://www.visioneng.com/products/software)

	ViCapture	DimensionOne	ViFox
Image capture	✓	✓	✓
Mark-up	✓	✓	✓
On-screen measurement		✓	✓
Reporting			✓
Advanced functions			✓

## For more information...

Vision Engineering has a network of offices and technical distributors around the world. For more information, please contact your Vision Engineering branch, local authorised distributor, or visit our website.

Distributor

Vision Engineering Ltd.  
(Manufacturing)  
Send Road, Send,  
Woking, Surrey, GU23 7ER, England  
Tel: +44 (0) 1483 248300  
Email: [generalinfo@visioneng.com](mailto:generalinfo@visioneng.com)

Vision Engineering Ltd.  
(Commercial)  
Monument House,  
Monument Way West,  
Woking, Surrey, GU21 5EN, England  
Tel: +44 (0) 1483 248300  
Email: [generalinfo@visioneng.com](mailto:generalinfo@visioneng.com)

Vision Engineering Inc.  
(Manufacturing & Commercial)  
570 Danbury Road,  
New Milford, CT 06776 USA  
Tel: +1 (860) 355 3776  
Email: [info@visioneng.com](mailto:info@visioneng.com)

Vision Engineering  
(Brasil)  
Email: [info@visioneng.com.br](mailto:info@visioneng.com.br)

Vision Engineering Ltd.  
(Central Europe)  
Anton-Pendele-Str. 3,  
82275 Emmering, Deutschland  
Tel: +49 (0) 8141 40167-0  
Email: [info@visioneng.de](mailto:info@visioneng.de)

Vision Engineering Ltd.  
(France)  
ZAC de la Tremblaie,  
Av. de la Tremblaie  
91220 Le Plessis Paté, France  
Tel: +33 (0) 160 76 60 00  
Email: [info@visioneng.fr](mailto:info@visioneng.fr)

Vision Engineering Ltd.  
(Italia)  
Via G. Paisiello 106  
20092 Cinisello Balsamo MI, Italia  
Tel: +39 02 6129 3518  
Email: [info@visioneng.it](mailto:info@visioneng.it)

Nippon Vision Engineering  
(Japan)  
272-2 Saedo-cho, Tsuduki-ku,  
Yokohama-shi, 224-0054, Japan  
Tel: +81 (0) 45 935 1117  
Email: [info@visioneng.jp](mailto:info@visioneng.jp)

Vision Engineering  
(China)  
Room 904B, Building B, No.970,  
Nanning Road, Xuhui Vanke Center  
Shanghai, 200235, P.R. China  
Tel: +86 (0) 21 5036 7556  
Email: [info@visioneng.com.cn](mailto:info@visioneng.com.cn)

Vision Engineering  
(South East Asia)  
P-03A-20, Impian Meridian,  
Jalan Subang 1,  
USJ 1, 47600 Subang Jaya,  
Selangor Darul Ehsan,  
Malaysia  
Tel: +604-619 2622  
Email: [info@visioneng.asia](mailto:info@visioneng.asia)

Vision Engineering  
(India)  
Email: [info@visioneng.co.in](mailto:info@visioneng.co.in)

Disclaimer – Vision Engineering Ltd. has a policy of continuous development and reserves the right to change or update, without notice, the design, materials or specification of any products. The information contained within this brochure/datasheet and to discontinue production or distribution of any of the products described.



Visit our multi-lingual website:

[www.visioneng.com](http://www.visioneng.com)